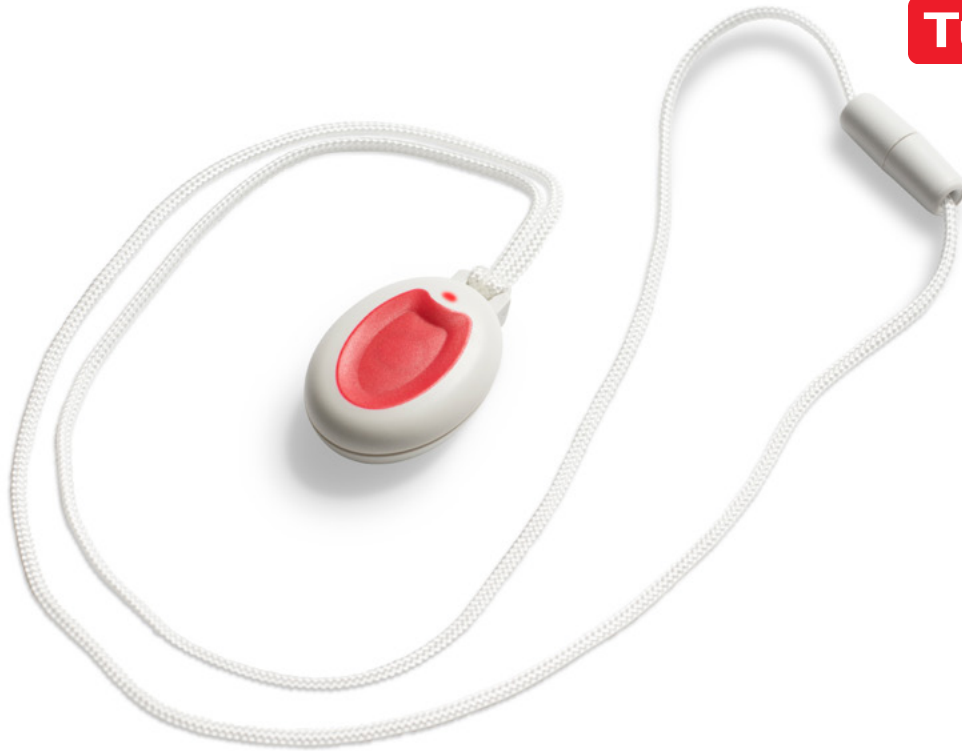




My Amie

What is it?

The MyAmie is a small, discreet, wearable pendant which allows the user to raise an alarm if they are in difficulty at home or in the garden.

Who is it for?

Suitable for anyone of any age living independently who may need to raise an alarm to request assistance, the MyAmie can offer peace of mind and an extra level of support. It can be particularly useful for older, infirm or disabled people, individuals who have recently been discharged from hospital or vulnerable people who could be at risk of crime or domestic abuse.

How does it work?

The waterproof pendant, which can be worn around the neck, on the wrist, or clipped to a belt, enables individuals to call for help by simply pressing a button. It can also be added to a key ring or mounted on the wall.

Within a 75m range of the Lifeline home unit, the MyAmie sends a radio signal to the Lifeline home unit or other Tunstall telecare enabled systems, which raises an alarm to the monitoring centre to take appropriate action in response.

Features

- Integrated alarm button with LED for visual reassurance that the alarm has been pressed
- A range of easy to wear options, including:
 - Wrist strap - with soft hypoallergenic, elasticated strap with clasp available in two widths
 - Neck cord - with safety break link
 - Belt clip - allowing the user to attach the pendant to their clothing
- Auto low battery monitoring and alerts

Benefits

- **Provides increased peace of mind and reassurance** for those living independently, their families and carers
- **Small, compact and comfortable design** is easy to wear and improves safety levels
- **Enables maximum reliability** operating on the 869MHz dedicated social alarm frequency

The MyAmie package

MyAmie is available as a single pendant on its own with a standard wrist band and neck cord. Optional extras are available separately. It is also included as part of the Smart Hub and the Lifeline Vi home solutions.

Why Tunstall?

We focus on using the latest digital and mobile technology to enable people to feel safe, secure and independent, giving them the freedom to live the life they choose. Our products combine secure digital connectivity and mobile platforms.

We help you provide...

- **Intelligent, unobtrusive, person-centred care.**
- **Personalised, proactive and predictive services to improve quality of life.**
- **Integrated health, housing and social care.**

For more information please visit:

uk.tunstall.com

Specification

Technical

Dimensions:

26 x 36 x 14mm (W x H x D)

Weight:

7g

Radio frequency:

European 869 MHz social alarm frequency

Radio range:

Up to 75m indoors

Up to 750m line of sight

Waterproof:

IP67 (Up to 1m for 30mins)

Compatibility:

All Tunstall 869MHz alarm systems

Standards

EMC:

EN 55032, EN 301 489-1, EN 301 489-3, EN 50130-4

Safety:

EN 60950-1

Radio:

EN 300 220-2

Social alarm:

EN 50134-2

Design, manufacture, installation and service: ISO 9001: 2008

Part numbers

MyAmie	P68005/01
Wrist Strap 20mm (pack of 10)	D6702145A
Wrist Strap 10mm (pack of 10)	D6702135A
Belt Clip (pack of 10)	D6702136A
Neck cord (pack of 10)	D6702137B
Easy press adapter	D6702112A



Tunstall Healthcare (UK) Ltd is a member of the Tunstall Group.

Contact us on: **t: 01977 661234** | **e: enquiries@tunstall.com** | **w: uk.tunstall.com** | **@TunstallHealth**

Our policy of continual development means that product specifications and appearance may change without notice. Tunstall does not accept any responsibility for any errors and omissions contained within this document. © 2018 Tunstall Group Ltd. Tunstall is a registered trademark. Tunstall Healthcare (UK) Ltd is a member of the Tunstall Group.

Version v10318